

ORAL HYGIENE



Create it 

ORAL HYGIENE

**NSK**

## Advanced oral care for dentists and hygienists

Removal of calculus and biofilm below and along the gum line together with tooth cleaning and polishing. Significantly improves the results of periodontal treatment.

NSK's wide range of oral hygiene products satisfies all these scenarios.

To maintain teeth in a healthy condition over a long period which is of great importance for patients.

NSK supports your daily clinical practice with effective and comfortable to use professional oral care equipment such as our high-performance ultrasonic and powder devices.

## INDEX

Varios Combi Pro .....	03 - 06
Ultrasonic .....	07 - 08, 13 - 14
Powder .....	21 - 22, 25 - 28
NSK Ultrasonic Solution .....	07 - 08
Varios Series .....	09 - 14
Ti-Max S970 .....	15 - 18
Prophy-Mate neo .....	19 - 20
FLASH pearl .....	21 - 22
Perio-Mate .....	23 - 28
Contra-angles .....	29 - 30
Spare Parts .....	31 - 32
Professional Care .....	33 - 34

# 2 in 1

Ultrasonic Scaling × Powder Therapy



For oral hygiene specialists.  
New from NSK, the global No.1 handpiece brand.

## Varios Combi Pro

## Two Maintenance Systems in One

Combining ultrasonic scaling and powder therapy for more effective and comfortable oral hygiene and biofilm control.



**Ultrasonic scaling for all**

**Powder therapy for supragingival**

**Powder therapy for subgingival (as option)**

**Innovative New Design Enables Clinicians to Safely and Efficiently Combine Two Applications.**

### Two systems in one

NSK has combined ultrasonic treatment and powder therapy into a single professional treatment unit that streamlines procedures and saves time and space.

### Automatically sets powder therapy functions for Supragingival or subgingival applications

The clinician simply attaches the relevant powder chamber, and the unit automatically selects the functions for supragingival or subgingival applications. This prevents selection errors.

### Independent unit and chamber

The detachable chamber and handpiece make it possible to easily switch between supragingival and subgingival applications and make the unit extremely user and maintenance friendly.

### No more clogging

By keeping the control unit, handpiece and handpiece hose separate the risk of powder blockages inside the unit has been completely eliminated. The components carrying powder and water are all easily accessible for maintenance and service.

### Keeps both powder types separate

Separate chambers, handpieces, and handpiece hoses for supragingival and subgingival procedures eliminate the risk of powder mixing or powder entering the control unit.


# Optimizing Manoeuvrability

## Clear Display

A clear function display enhances safety and efficiency by enabling clinicians to quickly confirm which programme is in use: supragingival, subgingival or ultrasonic, the 10 power levels and spray water volume.

## Large LCD Panel

The LCD panel offers excellent visibility and contrast and is specially coated to resist powder scratches.

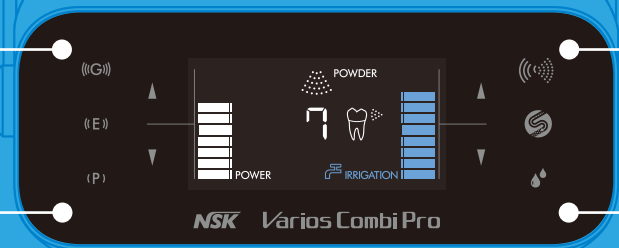


### During Ultrasonic function




### During Powder function



### When using Supragingival

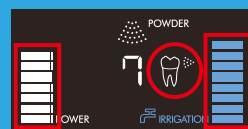

### When using Subgingival

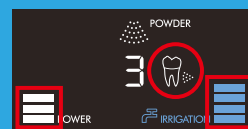

## Automatic Switching

An automatic recognition sensor switches from ultrasonic to powder function when a handpiece is picked up from the holder.

## Store Previous Settings

The system memory can store such treatment parameters as water volume, air pressure and ultrasonic power settings. Simply pick up the handpiece to switch settings.

\* For powder functions, settings change when replacing the chamber.

## Easy Maintenance

### Auto-cleaning function

Auto-cleaning automatically cleans the handpiece and hose after powder or ultrasonic applications. During the cleaning cycle, the display shows the time remaining to complete the process which takes approximately 45 seconds from start to finish.

### Easy maintenance

The nozzle, handpiece, handpiece hose and chamber are easy to disassemble.

### Smooth surfaces for easy cleaning

The body is seamless and smooth, making it easy to remove powder residues. The display material has been treated with a scratch resistant coating.

## Controls all Key Operations for Procedures

### Multi-Foot Control

NSK extensively researched foot movement and sensory feeling when developing the Multi-Foot Control to enable clinicians to operate the unit hands-free to facilitate uninterrupted procedures. Just using the foot the user can switch ultrasonic or powder functions on or off and change power settings. Combined with a hanger, this offers complete hands-free operations.



### Optional PERIO-SYSTEM

MODEL	ORDER CODE
Varios Combi Pro Perio Set	Y1003042

#### Contents

- PERIO-CHAMBER • PERIO-Hose • PERIO-Handpiece
- Perio Mate Nozzle Tip • Perio Mate nozzle tip remover • Perio Mate Powder

### Optional

Install the control unit and accessories on the cart. iCart Duo for VCP with the Water Bottle Set.

MODEL	ORDER CODE
iCart Duo	S9090V

- Dimensions: H 101.65 cm
- Weight: 16.5 kg



### Varios Combi Pro Basic Set

Optic	MODEL	ORDER CODE
●	Varios Combi Pro Basic Set	Y1004194

#### Contents

- PROPHY-Handpiece with 60° Nozzle • PROPHY-Hose • PROPHY-CHAMBER
- FLASH pearl (5 pcs. x 15 g sachets) • Maintenance kit • VA2LUX-Handpiece
- VA2LUX-Hose • 3 Scaling tips • Sterilization case
- Irrigation bottle (400 mL) • Foot control

#### Specifications

- Power supply: AC 100-240 V 50/60 Hz
- Frequency: 28-32 kHz
- Dimensions: W 225 x D 290 x H 162 mm (including chamber)
- Unit weight: 2.6 kg (Except attachment)
- Bottle volume: 400 mL






# iPiezo engine®

The unique ultrasonic AI engine responds intelligently to tooth surface conditions

"I'd like to be better at scaling", "I'd like an easier way to remove deposits"  
 We carefully considered the opinions of clinicians in the development of our unique next-generation AI ultrasound engine which is at the heart of the NSK ultrasonic series. The actions of the operator are transmitted accurately to the tip, as if it were an extension of the hand. This allows treatment to match tooth surface conditions. The result is faster and gentler treatment for patients and at the same time eco friendly as it operates efficiently with low power consumption.

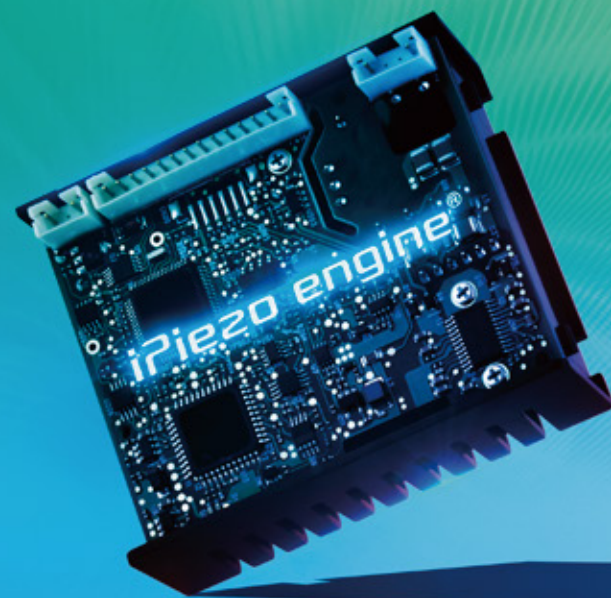

## Fitted with the next-generation ultrasonic AI engine, the iPiezo engine®

NSK's unique high-performance ultrasonic iPiezo engine allows for smooth operation as it transmits the operator's movements to the tip.

### Feedback Function


The NSK iPiezo engine supplies the appropriate power to the tip depending on the condition of the tooth surface. iPiezo engine provides a stable frequency output without dropping the power, even when removing hard calculus.

### Auto Tuning Function

The Auto Tuning Function provides stable frequency to the tip at the selected power range. It allows the operator to feel the power output which ensures the optimum result.


### Allows for steady scaling

Conventional ultrasonic scalers lose power when moving from soft to hard tartar. The load on the end of the tip is automatically detected by the action of the iPiezo engine allowing for steady scaling without power-loss.

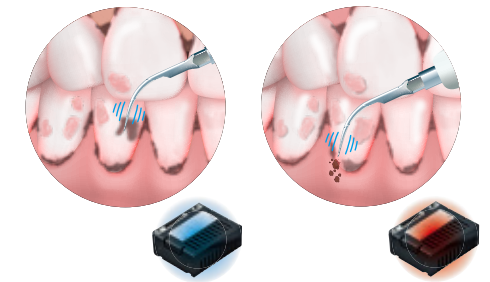

### LED A slim highly-ergonomic LED handpiece

The ø9.5 mm slim tip of the Varios2 LUX handpiece makes it easy to operate throughout the entire oral cavity. The slim handpiece facilitates a smooth approach to difficult to access areas such as the molars. This provides freedom of movement in the oral cavity allowing the tip to be applied at the correct angle. The easy to grip ultra-light handpiece reduces strain on the hands and fingers. With a wide spectrum of close to natural light LED illumination, the view is always bright with high recognition capability.

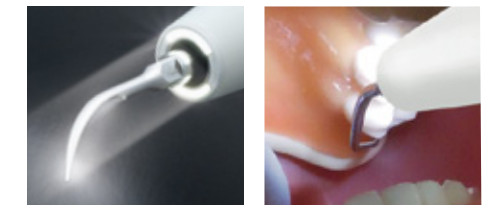

Autoclave sterilization and thermo-disinfection are both possible, delivering high safety standards




### Sterilization Case

The handpiece, tips and tip wrenches can be sterilized together for your convenience.



# Varios Series Be Smart.

NSK offers a wider range of progressively intelligent ultrasonic treatments with our unique next-generation AI engine and 70 varieties of tips

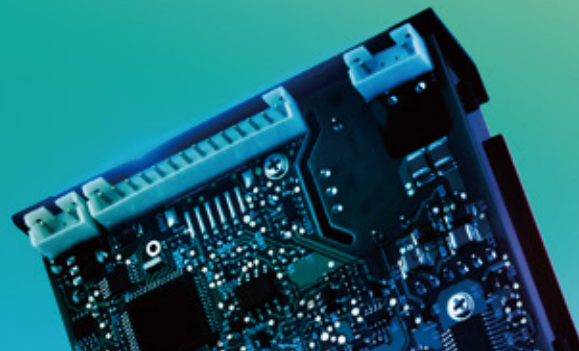

## Varios 970 Multifunction Portable Unit for Extreme Versatility

### The smart compact design provides high-performance

The functionality required for ultrasonic applications is available in the most-advanced model in our compact series, the Varios 970. Our high-capacity 2 bottle system allows users to focus on safe treatment without worrying about water irrigation. Simple daily maintenance with auto rinse function and seamless body design. Through intensive study of user feedback we have managed to achieve high operability and safety.

#### Large LED Display

The large LED display has a 10 stage numeric power display and bar graph to show the irrigation volume. Beyond this, intuitive layout of operating keys minimize user errors. Power adjustment and function selection can be performed easily and reliably.



#### Large Solution Bottles

The two bottles can be placed and removed easily and can be used for any type of irrigation to complement treatment. Large 400 mL bottles deliver ample irrigation capacity and reduced chair time. Repeated tests of our components and material studies ensure high durability even when used extensively with antiseptic irrigation solutions.



#### Improved Durability Against Antiseptic Solutions

The introduction of the resin based water pipe of the bottle, pump and handpiece, assures dramatically improved durability as it substantially reduces the corrosion and rust caused by the residuum of antiseptic solutions.




#### Auto Cleaning Mode

The Auto Cleaning Mode flushes the tubing, handpiece and tips at the push of a button when required. The two channels can be flushed independently.



#### Easy Replaceable Irrigation Pump

As you are able to replace the components of the irrigation system yourself, most issues can be dealt with on site.

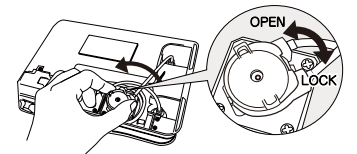

#### Varios 970 Complete Sets

Optic	MODEL	ORDER CODE
●	VA970 LUX	Y1001168
—	VA970	Y1001175

#### Contents

- Control unit • Handpiece • Handpiece cord • 3 Scaling tips (G4, G6, G8)
- Foot control (FC-70) • 2 x Irrigation bottles • 3 x Tip wrench with torque limiter
- AC power cord • Sterilization case

#### Specifications

- Power supply : AC 230 V 50/60 Hz
- Frequency : 28-32 kHz
- Max. output : 11 W (G mode)
- Dimensions : W 160 x D 270 x H 190 mm (Including bottle)
- Unit weight : 2.1 kg (Except attachment)
- Bottle volume : 400 mL (x2 bottles)



# Varios 370

## Compact Multipurpose Control Unit at Your Fingertips

Varios 370 has been especially developed to offer a compact, portable control unit allowing easy installation into any dental unit. It is simple to operate and suitable for a variety of clinical procedures.



### Advantages

- Powered by iPiezo engine
- Potable ultra compact body
- Simple and intuitive operation control panel

### Varios 370 Units Installation Examples



On the dental console arm with Stand (optional).

On/Under the dental bracket table with Control Unit Holder.

### Varios 370 Complete Sets

Optic	MODEL	ORDER CODE
●	VA370 LUX	Y1001388
—	VA370	Y1001391

**Contents**

- Control unit with handpiece cord • Control unit holder • AC adaptor
- Handpiece • Handpiece holder • Water filter set • Foot control (FC-71)
- 3 Scaling tips (G4, G6, G8) • Tip wrench with torque limiter

**Specifications**

- Power supply : AC 100-240 V 50/60 Hz
- Frequency : 28-32 kHz
- Max. output : 11 W (G mode)
- Dimensions : W 80 x D 115 x H 32 mm (without knob and cord)
- Unit weight : 0.43 kg (without attachment)



# Varios 170 Integrated Multi Function Ultrasonic Scaler

## The Smallest Module for Easy Integration into New and Existing Dental Equipment

Varios 170 is the world's smallest built-in optic control module. It features the NSK iPiezo engine which ensures stable power output all the times. The slimmer and extremely lightweight handpiece with twin LED is easily added to the existing handpiece delivery unit and is controlled via the existing foot control.



### Varios 170 Complete Sets

Optic	MODEL	ORDER CODE
●	VA170LUXS1	Y1001351
—	VA170S1	Y1001352

**Contents**

- Built in module • Handpiece • Handpiece cord • 3 Scaling tips (G4, G6, G8)
- 3 x Tip replacement wrench with torque limiter • Sterilization case


**Specifications**

- Power supply : AC 24 V ±10%
- Frequency : 28-32 kHz
- Max. output : 11 W
- Dimensions : W 63 x D 54 x H 26 mm
- Module weight : 59 g

# Tip Variation

## NSK Offers A Wide Choice of Ultrasonic Scaler Tips For Clinical Applications

There is a wide range of tip line up that can be selected by case and the state of the adhesive material.

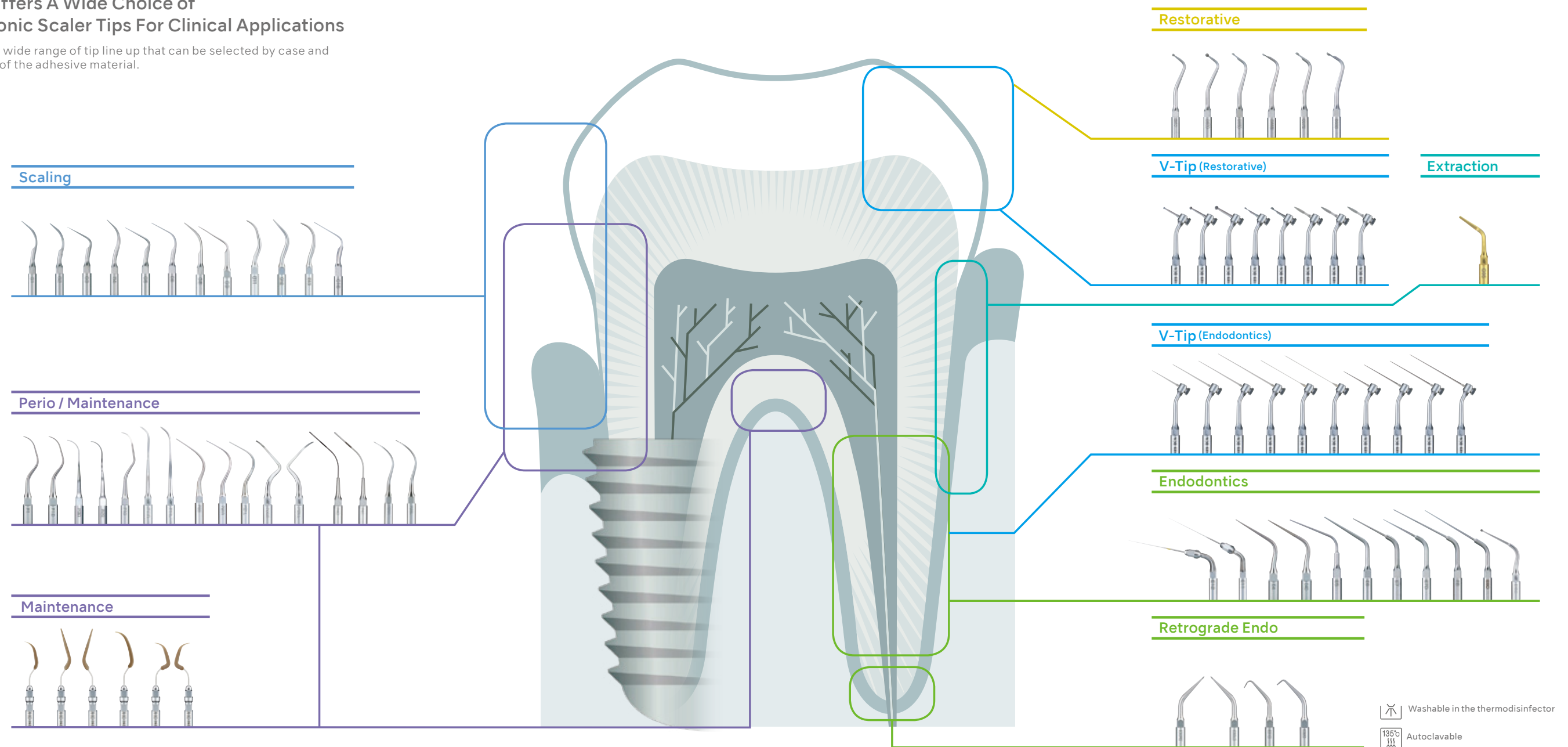

Washable in the thermodisinfectant  
Autoclavable

### Scaling / Perio tip cross section & characteristics of the tip types

Varios tips have distinct appearances, and also different cross-sections, component materials and surface finishes. The tip is selected by case and the state of the adhesive material. This makes for effective use without undue impact on the tooth surface or soft tissue.

Cross-section / tip	Material / tip end	Results of scaling treatment	Kind to teeth	For
G1, G5, G6 G8, G9, G16	Metallic mirror finish	⊙	△	Scaling at and below the gumline
P1, P20 P21R/L, P25R/L	Metallic mirror finish	△	⊙	Stain removal at and below the gumline / maintenance
P1D, P2D, P3D	Metal diamond machining	○	○	Scaling below the gum line (adjacent surface / bifurcation)

Cross-section / tip	Material / tip end	Results of scaling treatment	Kind to teeth	For
P10 P11R/L, P26R/L	Metallic sharp edge	⊙	△	Scaling below the gumline
P40, P41	Metallic, lateral tip, two point water irrigation	△	⊙	Irrigation
V-P10, V-P11R/L V-P12, V-P26R/L	Plastic	○	○	Prosthetic repair work Implants

See Tip Guide for details.





# Ti-Max S970 Air Scalers


## Air Scaler Evolution Supporting a Broad Spectrum of Care

The S970 air scaler offers a wide power range with easy swivel power adjustment and an extensive tip range for a wide variety of applications ranging from general scaling and root canal irrigation to minimally invasive procedures. The wide range of applications, comfortable grip and ease of maintenance makes the S970 air scaler the instrument of choice for dentists and hygienists alike. A titanium body, sleek and seamless design make the S970 really comfortable to hold even during long procedures whilst allowing for excellent visibility of the operating site at all times.




### Steady power supply to the tip

S970 air scalers provide stable power output and constant oscillating movement across the entire power range at approximately 6,000Hz. S970 air scalers are extremely gentle and therefore very popular in cases involving perio or high sensitivity.

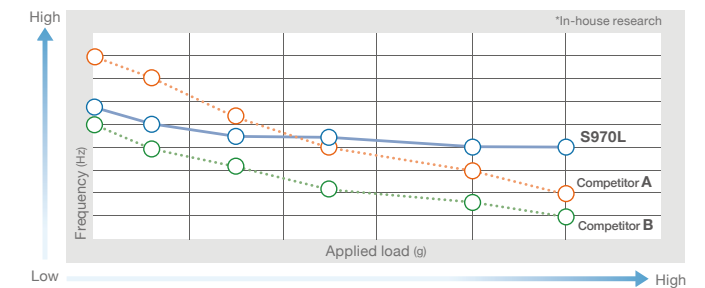

### Excellent visibility

The sleek design of the S970 allows for unrestricted visibility of the treatment area and excellent access to even the most hard to reach areas.

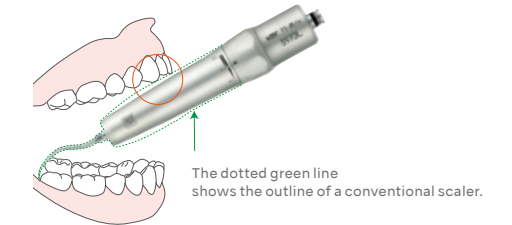

### Smooth and easy power control ring

NSK Power Ring offers three different power level settings for diverse clinical needs. While the gentle minimum mode is for delicate procedures, the maximum power mode is good for general scaling.






- Titanium Body with Scratch Resistant DURACOAT
- 3 level Power Ring • Cellular Glass Optics

Applications	Tips	Power Ring Setting		
		1	2	3
Scaling	S1, S2, S3	[Full range arrow]		
Perio	S10, S20, S21R, S21L	[Low to High arrow]		
Perio (Maintenance Tip)	S35+V-Tip (V-P10, V-P11R, V-P11L, V-P12)	[Low to High arrow]		
Restorative (For Preparation of crown cores)	S81D, S82D, S83D, S86D, S71D, S72D	[Low to High arrow]		
Restorative (For Minimal Intervention)	S65D, S66D, S67D, S68D	[Low to High arrow]		
V-TIP SYSTEM (Endodontics)	S75+V-Tip (V-U15-V-U35)	[Low to High arrow]		

The above power level is standard. Please start use from a low level to correspond to the symptoms.

### Coupling Lineup

	NSK		KaVo®	Dentsply Sirona®
MODEL	<b>S970L</b>	<b>S970</b>	<b>S970KL</b>	<b>S970SL</b>
ORDER CODE	<b>T1005</b>	<b>T1020</b>	<b>T1013</b>	<b>T1015</b>
Optic	●	—	●	●
COUPLINGS				
	NSK Coupling		KaVo® Multiflex®	Sirona® Quick Coupling

- Contents**  
 • 3 tips (S1, S2 & S3) • Tip wrench (CR-20) • Tip cover (L)

 Washable in the thermodisinfectant  Autoclavable

### AS2000 Air Scalers

- Frequency of oscillation : 5,600-6,000 Hz • 3 tips (S1, S2 & S3)
- Tip Wrench (CR-20) • Tip cover (L)




Optic	MODEL	ORDER CODE
—	<b>AS2000 M4</b>	<b>T698</b>

• For Midwest 4 hole

# Tip Variation

## NSK Offers A Wide Choice of Air Scaler Tips For Clinical Applications


There is a wide range of tip line up that can be selected by case and the state of the adhesive material.

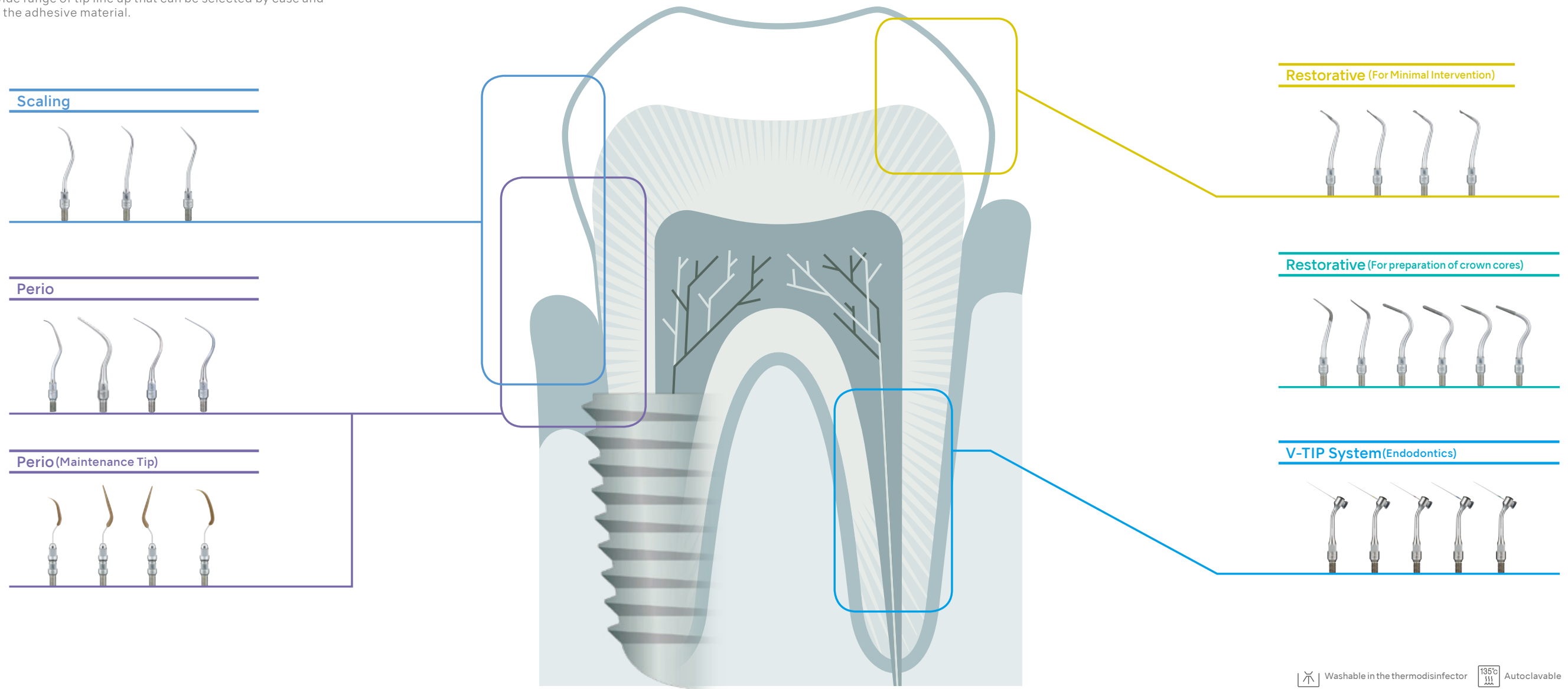

Washable in the thermodisinfectors 135°C Autoclavable

TI-Max S970

### Scaling

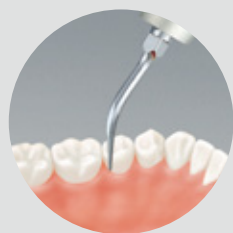


Removal of sub gingival calculus

### Perio




Removal of gingival biofilm

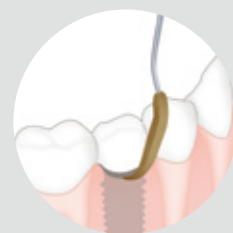


Maintenance of prostheses and implants

### Detailing and finishing

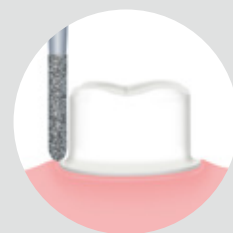

Used for the final finish margins

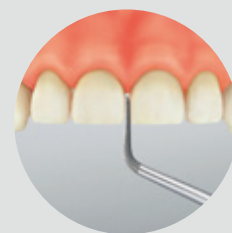

Preparation of interproximal areas

The tips can be used for finishing margins after prepping the area with a conventional dental bur. That enables you to finish minute sharp edges and gingival sulcus without damage, something that cannot be accomplished with rotary instruments alone. The torpedo shaped tip is suitable for the final formation of an adjacent surface margin such as for laminated veneers. Low risk of excessive grinding and dentin exposure enables more conservative formation of the adjacent surfaces.

See Tip guide for details.





# Prophy-Mate neo Air Powered Tooth Polishing System

## Prophy-Mate neo, the best polishing tool for beautiful white teeth

Supragingival powder polishing is highly effective when used in combination with "FLASH pearl." NSK's unique twin nozzle achieves powerful powder discharge. The long and slim handpiece in conjunction with the small nozzle tip dramatically improves accessibility around all tooth surface. Cleaning can be performed quickly even in interproximal areas teeth and in usually difficult to reach areas. The ergonomically designed grip is comfortable to hold and allows for flexible operation.



### Improved Polishing Power

With NSK original twin nozzles, greater polishing power is realized. With powerful jet, cleaning time is greatly shortened thus reduces treatment time.



### Easy Detachable Nozzles with 60° and 80°

Selectable 60° and 80° nozzles maximize the cleaning ability. 60° nozzle works especially well for molar area and occluding surface and 80° nozzle for anterior tooth.



### Smooth, Easy 360° Swivel

Twin joints rotate smoothly 360 degrees for easy operability, responding flexibly to minute fingertip movements with effortless adjustment of the nozzle angle.










### Simple One-touch Connection

The Prophy-Mate neo is designed for a one-touch attachment to NSK couplings and other major coupling brands. The handpiece also can be attached to the powder case with one-touch operation.



### Coupling Lineup

	NSK	KaVo®	Dentsply Sirona®	W&H®	Bien Air®	M4	B2
MODEL	PMNG-PTL-P	PMNG-KV-P	PMNG-SR-P	PMNG-WH-P	PMNG-BA-P	PMNG-M4-P	PMNG-B2-P
ORDER CODE	Y135029	Y135030	Y135031	Y135032	Y135033	Y135034	Y135035
COUPLINGS	 NSK Coupling	 KaVo® Multiflex®	 Sirona® Quick Coupling	 W&H® Roto Quick®	 Bien Air® UNIFIX®	 Connect direct to the handpiece tubing MIDWEST	 Connect direct to the handpiece tubing BORDEN

**Contents**  
 • Handpiece with 60° nozzle head • Powder case • Replacement nozzle head 80°  
 • 5 pcs. x 15 g sachets of FLASH pearl cleaning powder • Maintenance kit

\*No light during use (when using optic coupling).

 The handpiece can be washed in the thermosymbol

 Autoclavable



# FLASH pearl Cleaning Powder

## Outstanding Cleaning Results with Gentle Powder


Not salty! The powder has less impact due to the delicate flavour and provides patients with gentle tooth surface cleaning. The secret behind intense cleaning and a smooth finish? Calcium-based spherical particles. Fine powder with an average particle size of 54µm move as if rolling across the surface of the tooth, limiting damage to the tooth surface and rapidly removing plaque and staining. The "FLASH pearl" treatment is effective in a short space of time which encourages patients to schedule regular cleaning appointments.



## Sparkling results every time

FLASH pearl, a cleaning powder for Prophy-Mate neo and Varios Combi Pro, is a spherical particulate powder designed to roll smoothly over tooth surfaces and penetrate all corners for rapid cleaning. The spherical shape of the powder granules lessens the likelihood of damage to teeth and soft tissue while removing stains and plaque quickly and safely.


☆ Adverse reactions to this product are rare.

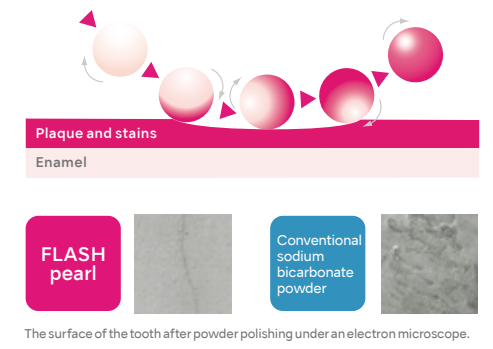

## Fast cleaning in tight spaces

The nozzle is held close to the tooth surface at an angle close to horizontal (10° to 60°) and maneuvered in an intricate motion as if the powder is moving across the tooth surface. The wide discharge angle and freer range of discharge motion allow for smooth operation and contribute to improved treatment efficiency.

☆ Please do not apply directly to the gumline.

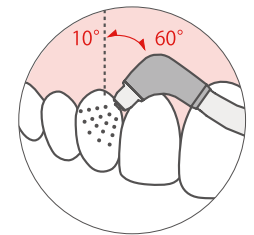

## Zero salinity, gentle ingredients with low impact on patients

As the powder is composed of 94% calcium carbonate, it is not unpleasant for your patients and absent in any salt-like taste. There's no need to worry about patients on restricted salt diets.

## Clogging is also minimized

The fine powder particles don't dissolve easily in water so clogging is reduced. Maintenance of equipment is easy as there is less powder build-up inside instruments.

MODEL	ORDER CODE
FLASH pearl (Bottle)	Y900698

- 4 bottles packed in a box
- Comes with a measuring cup for your convenience
- Each bottle contains 300 g



# Perio-Mate Biofilm Eraser

## Extensive range of applications ranging from supragingival to subgingival use

By removing the Perio Mate Nozzle Tip, the subgingival area up to 3 mm from the gum line that requires the most frequent maintenance can be treated. In addition, there is a powder flow rate adjustment function to match operational site conditions, to enable a wide range of possible treatments from fast precise removal of biofilm to gentle and efficient biofilm removal around implants and in the delicate periodontal pocket.



Powder Pro Care, Above And Below The Gums

Efficient professional care that's kind to teeth. It is easy to install powder therapy alongside existing facilities.

- The "handy type" can be connected to the chair unit during use in one action.
- The nozzle tip is designed for safety and operability.
- The powder flow channel is short and direct with low risk of blockage.
- Fast and effective biofilm removal.

### Powder flow volume adjustment

The Perio-Mate's convenient Powder Control Ring makes fine-tuning powder flow simple and intuitive with quick fingertip adjustments. In addition to easily controlling and stopping the flow of powder in this way, you can cleanse the treatment area from biofilm uninterruptedly.



### Smooth 360-degree rotation for safe and accurate operation

The handpiece can easily be rotated using the fingertips even while performing complex operations. The nozzle tip allows easy access to difficult to reach areas, and does not adversely affect neighbouring teeth and contours.



### Slim nozzle facilitates operation in the narrow confines of the oral cavity

The long, narrow nozzle shape is designed to allow easy access to areas that require inserting the handpiece at an angle, such as between molars.










### Product maintenance

The user can easily disassemble the nozzle, handpiece and powder case. Powder clogging is unlikely when handled and maintained correctly but is easy to manage by following the maintenance guidelines and the relevant accessories included with the product. All the parts of the handpiece (other than the nozzle tip) can be repeatedly cleaned ultrasonically or sterilized in an autoclave.



### Coupling Lineup

	NSK	KaVo®	Dentsply Sirona®	W&H®	Bien Air®	M4	B2
MODEL	PER-PTL-P	PER-KV-P	PER-SR-P	PER-WH-P	PER-BA-P	PER-M4-P	PER-B2-P
ORDER CODE	Y1002654	Y1002655	Y1002656	Y1002657	Y1002658	Y1002659	Y1002660
COUPLINGS	 NSK Coupling	 KaVo® Multiflex®	 Sirona® Quick Coupling	 W&H® Roto Quick®	 Bien Air® UNIFIX®	 Connect direct to the handpiece tubing MIDWEST	 Connect direct to the handpiece tubing BORDEN

Contents  
 • Handpiece • Powder case • Nozzle tip (40 pcs) • Maintenance kit • Perio Mate Powder 160 g x1

\*No light during use (when using optic coupling).

 The handpiece can be washed in the thermodisinfecteur


 Autoclavable



## Supragingival

Removal of plaque biofilm from the gingival margin (up to approximately 3 mm inside the pocket) \*Detach nozzle tip


### Removal of plaque biofilm at interdental and gingival margin

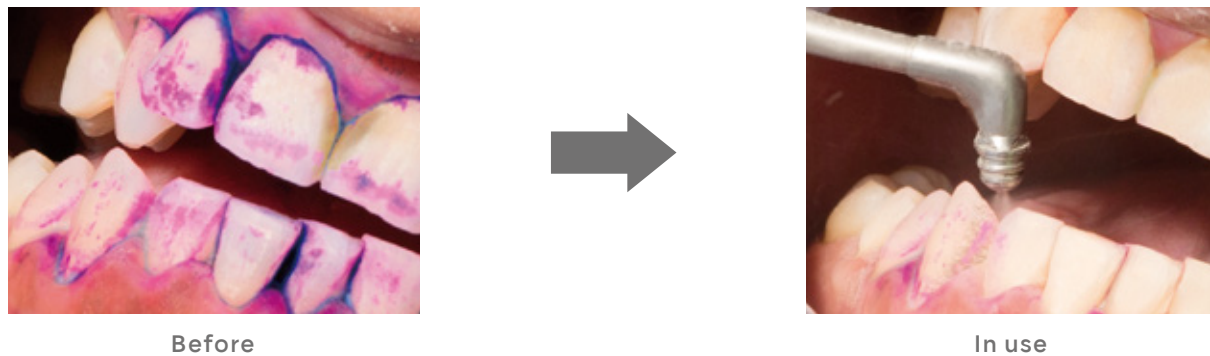

- Easier to remove biofilm plaque and minor stain around supragingival and gingival margin.
- Enable to remove plaque and biofilm from crowded dentition, interdental and pit and fissure, where hardly accessed with rubber cups.
- Removal of minor stain attached extensively.
- Cleaning around dental implants and prosthetic restorations.
- Less of burden to teeth and gums due to non-contact method.
- Less of burden to prosthetic restorations.

### Fast and effective removal capability

Biofilm on teeth surface can be removed in app. 50.5 seconds per 1/4 dentition.

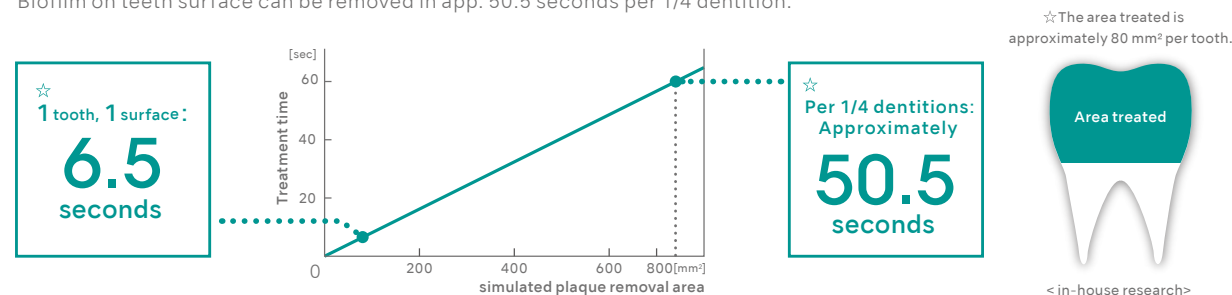


### Observation of the treated area at 200-fold with a confocal laser microscope\*


Surface of prosthetic restorations can be treated without damages.

\*CLSM Lasertec


#### Hybrid hard resin

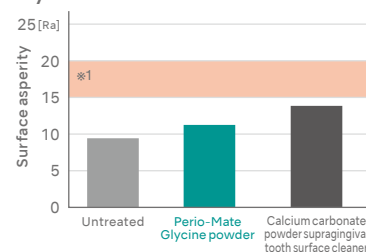

#### Gold and silver palladium alloy

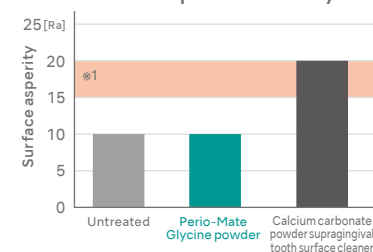


#### Titanium alloy

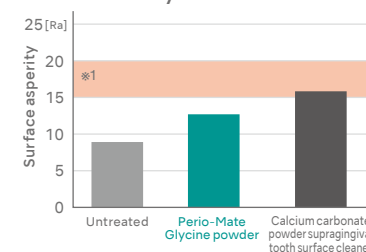

\*1 Average surface asperity of regular diamond point


Source: Oral Health Engineering, School of Oral Health Care Sciences, Tokyo Medical and Dental University

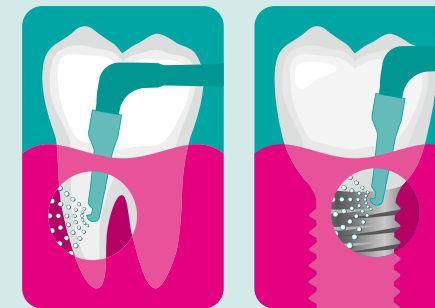

## Subgingival

Erasing biofilm in periodontal pocket from 3 - 6 mm below the gingival margin


### Removal of biofilm in dental pocket

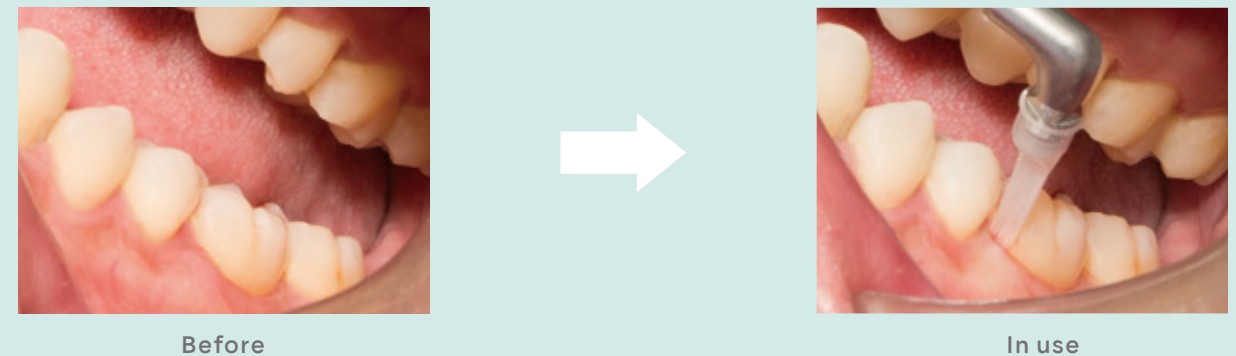

- Biofilm in the periodontal pocket can be removed in approximately 60 seconds per quadrant.
- Deposits can be removed without direct contact with almost none of the extensive biofilm and plaque adhesion remaining and the results are effective. It is possible to clean areas the tip cannot reach with the powder injection effect.
- Optimized the flow of powder based on fluid analysis.
- Compared to other maintenance equipment, less biofilm and plaque remains with wide powder dispersal, even in areas that cannot be seen.
- Smooth powder injection ensures gentle treatment in the periodontal pocket.
- Flexible plastic nozzle tip design ensures accessibility and is gentle to the pocket.

### Fast and effective removal capability

Biofilm in the periodontal pocket can be removed in approximately 60 seconds per quadrant.

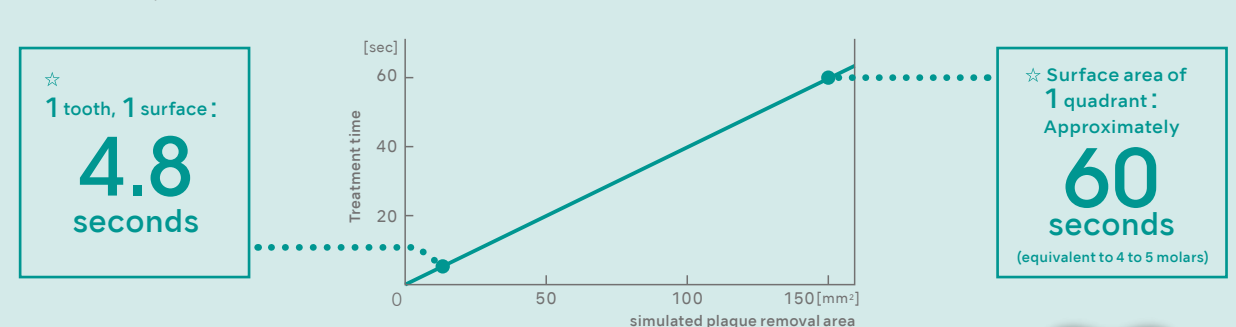

1. Used for SPT (supportive periodontal therapy) for pockets around 6mm that result from periodontal disease.\*2
2. Spraying each tooth for 20 seconds (5 seconds per surface) is best to eradicate subgingival biofilm inside delicate periodontal pockets.\*2
3. Reciprocal motion of the nozzle tip and fine circular movements are vital to remove biofilm effectively.



\*2. Flemmig TF, Hetzel M, Topoll H, Gerss J, Haeberlein I, Petersilka G., Subgingival debridement efficacy of glycine powder air polishing., J Periodontol. 2007 Jun;78(6):1002-10., Petersilka GJ, Steinmann D, Häberlein I, Heinecke A, Flemmig TF., Subgingival plaque removal in buccal and lingual sites using a novel low abrasive air-polishing powder., J Clin Periodontol. 2003 Apr;30(4):328-33.




# Performance & Technology

## Reduction of bacteria


The air-polish process with Perio Mate Powder (glycine-based) results in significantly better bacterial reduction than manual plaque removal.

Petersilka GJ, Tunkel J, Barakos K, Heinecke A, Haberlein I, Flemmig TF. Subgingival plaque removal at interdental sites using a low-abrasive air polishing powder. J Periodontol. 2003 Mar;74(3):307-11.

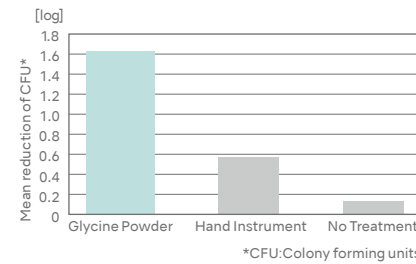


## Powder particle comparisons

The fine Perio Mate Powder is ideal for periodontal pocket treatment.

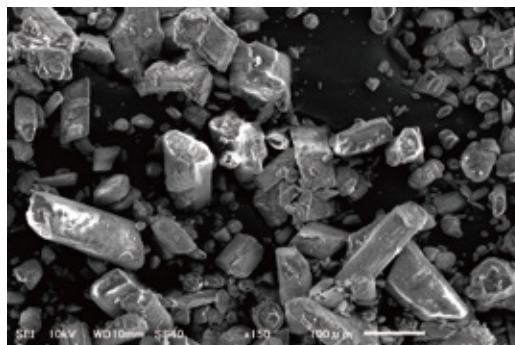


**Conventional sodium bicarbonate powder**

Average particle diameter: 73 µm  
Scanning electron microscope image magnification of 150x

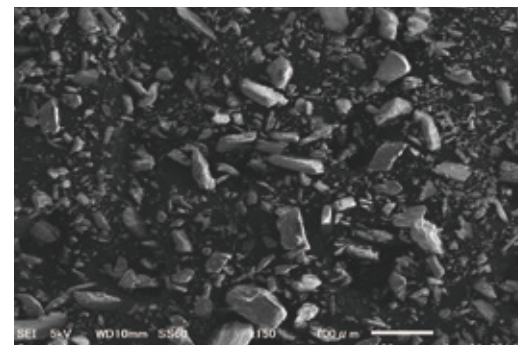

**Perio Mate Powder (glycine)**

Average particle diameter: 25 µm  
Scanning electron microscope image magnification of 150x

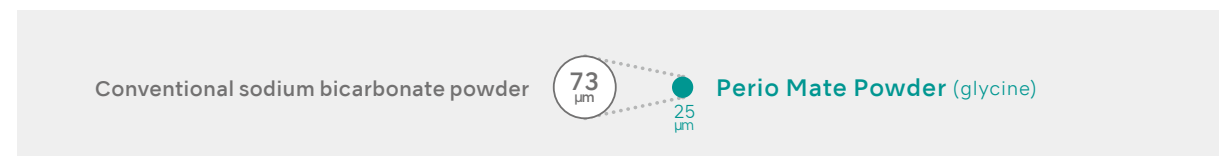

## Perio Mate Powder

- Subgingival airpolishing powder.
- Principal ingredient is glycine, an amino acid present in proteins.
- Highly water soluble to prevent subgingival retention.
- Lower Mohs hardness than dentin. \*Glycine ranks 2 on the Mohs hardness scale, while dentin is 2 to 2.5.
- Average powder particle diameter is 25 µm.

MODEL	ORDER CODE
<b>Perio Mate Powder</b>	<b>Y900938</b>

**Contents**  
• Two bottles (160 g each)



## NSK's innovative nozzle tip

NSK's new nozzle tip design incorporates semitransparent materials that promote greater visibility during procedures. Hygienic single use tips mean enhanced usability and safety.




## Flexible plastic nozzle tip design ensures accessibility and is gentle to the pocket

The slim nozzle tip possesses just the right degree of flexibility to pass smoothly over tooth surfaces, and can be easily inserted into tight pockets and the root surfaces of adjacent teeth without causing the patient discomfort.



## Optimized using powder flow analysis

Turbine development technology uses accumulated fluid analysis which allows 'spraying' that takes into account the powder flow in the pocket. It is not mono-directional, rather the flow is even and all-directional, minimising soft-tissue damage and subgingival plaque and biofilm can be removed efficiently. With respect to potential damage, powder (air) and water are separately discharged onto periodontal tissue.

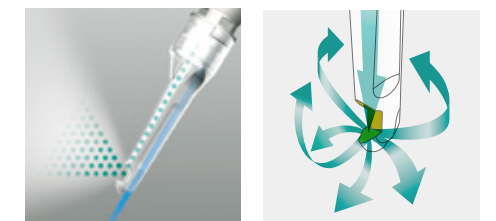

## Perio Mate Nozzle Tip

- Disposable
- Sterile
- Individually packaged



MODEL	ORDER CODE
<b>Perio Mate Nozzle Tip</b>	<b>Y1002741</b>

**Contents**  
• 40 nozzles



# Contra-angles Prophylaxis Handpiece



## Screw cup polishing

### Ti-Max X



Optic	MODEL	ORDER CODE
☉	<b>X57L</b>	<b>C609</b>

- 16:1 Reduction • Max Speed : 2,500 min<sup>-1</sup>
- Titanium body & DURACOAT • Screw-in cups and brushes
- 360° rotation



Special dust proof

### FX series



Optic	MODEL	ORDER CODE
—	<b>FX57m</b>	<b>C1054001</b>

- 4:1 Reduction • Max Speed : 5,000 min<sup>-1</sup> • Stainless steel body
- Screw-in cups and brushes • 360° rotation



Paste defence system

### EX series



#### Heads

MODEL	ORDER CODE
<b>AR-Y (S)</b>	<b>C187</b>

- Max Speed : 5,000 min<sup>-1</sup>
- For screw-in cups and brushes
- 360° rotation

#### Shanks

MODEL	ORDER CODE
<b>EC-30M</b>	<b>C100007</b>

- 1:1 Direct Drive

MODEL	ORDER CODE
<b>AR-Y (K)</b>	<b>C188</b>

- Max Speed : 5,000 min<sup>-1</sup>
- For snap-on cups • 360° rotation

MODEL	ORDER CODE
<b>ER4M</b>	<b>C634002</b>

- 4:1 Reduction

## EVA tip polishing

### Ti-Max X



Optic	MODEL	ORDER CODE
☉	<b>X55L</b>	<b>C608</b>
—	<b>X55</b>	<b>C618</b>

- 4:1 Reduction • Max Speed : 10,000 min<sup>-1</sup>
- Titanium body & DURACOAT
- For Profin Tip / EVA Tip / Orthodontics File (Shank ø2.36-2.40 mm)\*
- Vertical Reciprocating (Stroke Length 1.4 mm) • Single Spray
- \*Metallic shank cannot be used.



With tip lock mechanism

### EX series



#### Head

MODEL	ORDER CODE
<b>EVA-Y</b>	<b>C235</b>

- Max Speed : 10,000 min<sup>-1</sup> • Eva tips
- Vertical reciprocating (stroke length 1.4 mm)

#### Shank

MODEL	ORDER CODE
<b>ER4M</b>	<b>C634002</b>

- 4:1 Reduction

# Spare Parts

## Varios Combi Pro

MODEL	PRODUCTS	ORDER CODE
<b>Chamber</b> For Prophy use		U1135018
<b>Chamber</b> For Perio use		U1135019
<b>Handpiece Hose</b> For Prophy use		T1051001
<b>Handpiece Hose</b> For Perio use		T1053001
<b>Powder Handpiece</b> For Prophy use		T1074001
<b>Powder Handpiece</b> For Perio use		T1077001
<b>Nozzle</b> For Prophy use		T1074012
<b>Nozzle</b> For Perio use		T1077012
<b>VCP Bottle Set</b> Bottle and Cap		U1135017
<b>VA2-LUX-HP</b> Optic Handpiece		E351050
<b>VA2LUX-SC</b> Handpiece Cord for VA2LUX		E351054

MODEL	PRODUCTS	ORDER CODE
<b>FC-79</b> Foot Control		Z1256001
<b>Water Filter Set</b> 3.6 m		U1135078
<b>Water Pressure Reducing Valve</b>		U1135093
<b>Irrigation Pump</b>		U1135079
<b>O-ring set</b> Chamber Base Connector (Pack of 5)		Y1003435
<b>Chamber Cap</b>		U1135089
<b>O-ring</b> For Chamber Cap		Y1003545
<b>O-ring set</b> Irrigation Bottle		Y1003544
<b>O-ring</b> 1 for the Handpiece Cord		Y1003546
<b>O-ring set</b> Handpiece Hose		Y1003500
<b>O-ring set</b> Nozzle		Y1003501

MODEL	PRODUCTS	ORDER CODE
<b>Air Hose</b> 3.6 m		D0914406 036
<b>Ring Wrench</b>		T1074601
<b>Perio Mate Nozzle Tip</b> (Pack of 40)		Y1002741
<b>Nozzle Tip Remover</b>		T1043751
<b>Lubrication Oil</b>		Z016112
<b>Cleaning Wire (L)</b>		T798091
<b>Cleaning Wire (S)</b>		T798095
<b>Cleaning File</b>		T798092
<b>Cleaning Brush</b>		T798093
<b>Unlocking tool</b>		U1135256
<b>Fuse</b>		D1202220 010

## Varios 970 / 370 / 170

MODEL	PRODUCTS	ORDER CODE
<b>VA2-LUX-HP</b> Optic Handpiece		E351050
<b>VA2-HP</b> Non-Optic Handpiece		E350050
<b>VA970LUX-SC</b> Handpiece Cord for Varios 970 LUX		E351055
<b>VA970-SC</b> Handpiece Cord for Varios 970		E350055
<b>VA170LUX-SC</b> Handpiece Cord for Varios 170 LUX		E351056
<b>VA170-SC</b> Handpiece Cord for Varios 170		E350056
<b>FC-70</b> Foot Control (2.5 m)		Z1008001
<b>FC-71</b> Foot Control (5 m) For Varios 370		Z1008002
<b>AC Adaptor</b> (5 m) For Varios 370		Y1001792

MODEL	PRODUCTS	ORDER CODE
<b>Water Supply Connector</b>		U387030
<b>Water Filter Set (2 m)</b> 2 water supply tube, water filter, water filter case		U387040
<b>Water Filter Set (4 m)</b> For Varios 370 2 water supply tube, water filter, water filter case		U1007080
<b>Water Supply Tube (2 m)</b>		U387580
<b>Water Supply Tube (4 m)</b> For Varios 370		D0914406 040
<b>Water Filter Case</b>		Y900069
<b>Tube Ring</b> (Pack of 2)		Y900077
<b>Handpiece Holder</b> For Varios 370		10001505
<b>VA Bottle 400</b> (400 mL) For Varios 970		Z1047002

MODEL	PRODUCTS	ORDER CODE
<b>Irrigation Pump</b> For Varios 970		10000643
<b>O-ring</b> For Handpiece Cord		D0310020 080
<b>O-ring</b> Thin section For VA Bottle 400		D0312090 100
<b>O-ring</b> Thick section For VA Bottle 400		D0310075 150
<b>Varios 350/370 Stand</b>		Z254
<b>Varios 350/370 Mounting Bracket</b>		Z259
<b>Control Unit Holder</b> For Varios 370		Y1001812

## Varios Combi Pro / Varios 970 / 370 / 170

MODEL	PRODUCTS	ORDER CODE
<b>Water Filter</b>		U387042
<b>CR-10</b> Tip Wrench		Z221076
<b>E Tip Replacement Wrench</b>		Z217399

MODEL	PRODUCTS	ORDER CODE
<b>Tip Holder</b>		Z221A080
<b>Tip Cover L</b>		Z217850
<b>Tip Cover S</b>		Z217851

MODEL	PRODUCTS	ORDER CODE
<b>Sterilization case</b>		Z1035001
<b>Spanner Wrench</b> 5x8 (Pack of 2)		Y1001301

## Prophy-Mate neo

MODEL	PRODUCTS	ORDER CODE
<b>PMN-HP</b> Handpiece with 60° nozzle head		T797010
<b>Nozzle 60°</b>		T797012
<b>Nozzle 80°</b>		T797013
<b>O-ring</b> neo Rear Nozzle		D0312029 070

MODEL	PRODUCTS	ORDER CODE
<b>Blower Nozzle-PTL</b> neo		T798094
<b>Blower Nozzle-KV</b> neo		T886094
<b>Blower Nozzle-SR</b> neo		T887094
<b>Blower Nozzle-WH</b> neo		T888094

MODEL	PRODUCTS	ORDER CODE
<b>Blower Nozzle-BA</b> neo		T889094
<b>Blower Nozzle-M4</b> neo		T884094
<b>Blower Nozzle-B2</b> neo		T885094

## Perio-Mate

MODEL	PRODUCTS	ORDER CODE
<b>PER-HP</b>		T1043010
<b>Nozzle</b>		T1043012
<b>O-ring</b> PER Rear Nozzle		D0312023 060
<b>Nozzle Tip Remover</b>		T1043751

MODEL	PRODUCTS	ORDER CODE
<b>Blower Nozzle-PTL</b> PER		T1042094
<b>Blower Nozzle-KV</b> PER		T1044094
<b>Blower Nozzle-SR</b> PER		T1046094
<b>Blower Nozzle-WH</b> PER		T1048094

MODEL	PRODUCTS	ORDER CODE
<b>Blower Nozzle-BA</b> PER		T1050094
<b>Blower Nozzle-M4</b> PER		T1052094
<b>Blower Nozzle-B2</b> PER		T1054094

## Prophy-Mate neo / Perio-Mate

MODEL	PRODUCTS	ORDER CODE
<b>Powder Case Cover</b>		T798133
<b>Gasket</b>		T798135
<b>O-ring</b> Front Nozzle		D0312039 070
<b>Ring Wrench</b>		T797601

MODEL	PRODUCTS	ORDER CODE
<b>PTL O-ring Set</b> (Pack of 5)		Y900580
<b>PM O-ring Set</b> (Pack of 3)		Y900054
<b>Lubrication Oil</b>		Z016112
<b>Cleaning Wire (L)</b>		T798091

MODEL	PRODUCTS	ORDER CODE
<b>Cleaning Wire (S)</b>		T798095
<b>Cleaning File</b>		T798092
<b>Cleaning Brush</b>		T798093

## Ti-Max S970 / AS2000 Air Scaler

MODEL	PRODUCTS	ORDER CODE
<b>CR-20</b> Tip Replacement Wrench with Torque Limiter		Z1090001
<b>E Tip Replacement Wrench</b>		Z217399

MODEL	PRODUCTS	ORDER CODE
<b>Tip Holder</b>		Z252A080
<b>Tip Cover L</b>		Z217850

MODEL	PRODUCTS	ORDER CODE
<b>Tip Cover S</b>		Z217851
<b>S970 O-ring Set</b> For S970 (Pack of 4)		Y1001953



# Professional Care

## Selecting the right equipment is the first step in professional care

Professional care involves a variety of treatments such as the prevention and improvement of periodontal disease, subgingival scaling and the removal of pathogenic bacteria, the removal of biofilm and prevention of dental caries through prevention of plaque re-adhesion and the removal of calculus and stains to restore the teeth to their original brilliance. In order to provide total professional care, it is important to select quality equipment to provide a suitable treatment and operational efficiency to suit each patient's oral care requirements.



KaVo and Multiflex are registered trademarks of KaVo Dental GmbH.  
Dentsply Sirona and Sirona are registered trademarks of Dentsply Sirona Inc.  
W&H and Roto Quick are registered trademarks of W&H Deutschland GmbH.  
Blen Air and UNIFIX are registered trademarks of Blen-Air Holding SA.

**NAKANISHI INC.** [www.nsk-dental.com](http://www.nsk-dental.com)  
700 Shimohinata, Kanuma, Tochigi 322-8666, Japan

**NSK Oceania Pty. Ltd.** [www.australia.nsk-dental.com](http://www.australia.nsk-dental.com)  
Unit 12, 809-821 Botany Road, Rosebery NSW 2018, Australia

**NSK Oceania Ltd.** [www.new-zealand.nsk-dental.com](http://www.new-zealand.nsk-dental.com)  
Ground Floor, Suite D, 272 Parnell Road, Parnell, Auckland 1052, New Zealand

Specifications are subject to change without notice.